


AUXILOCK[®] GYROLOCK PEEK OPTIMA KNOTLESS SCREW IN ANCHOR

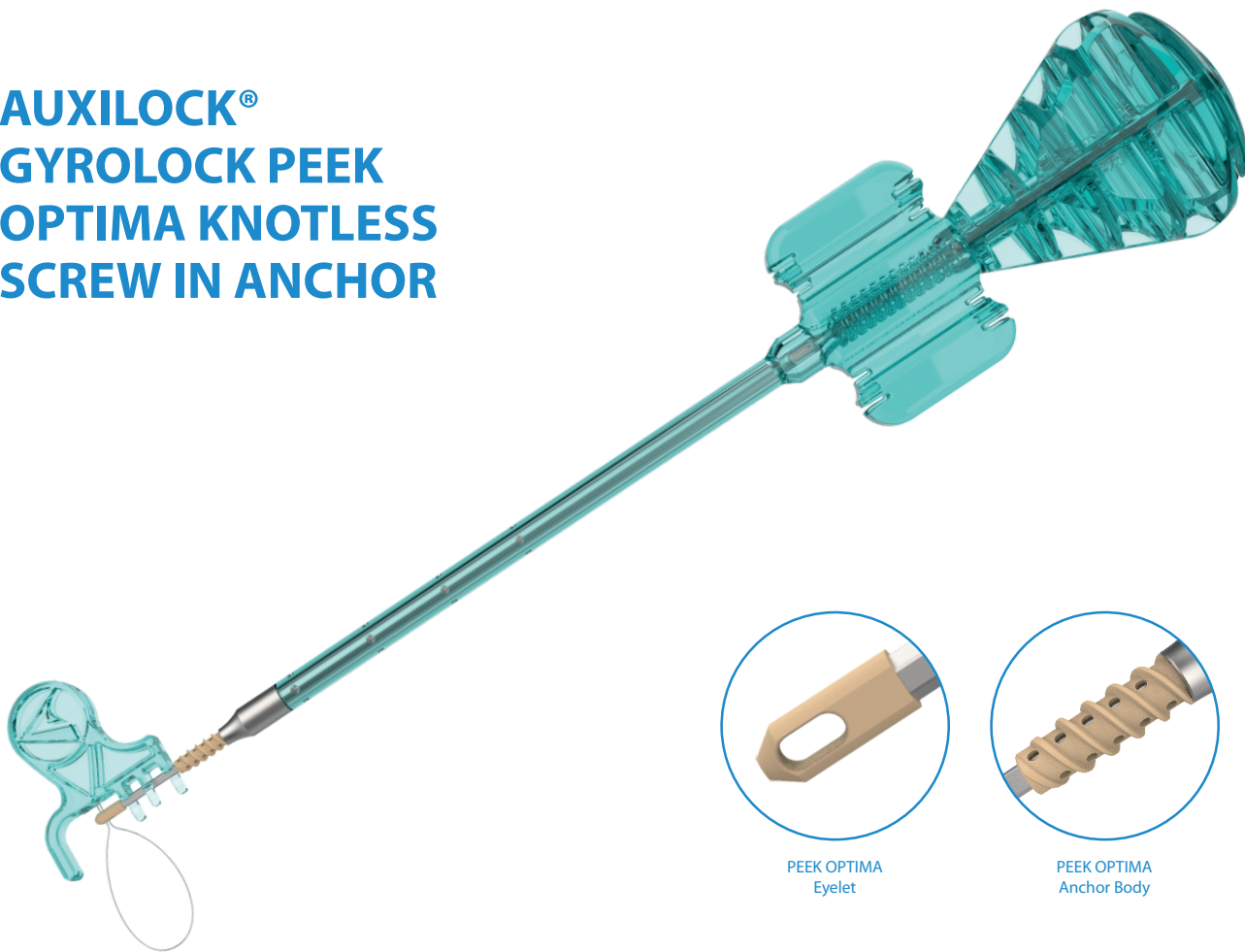

GyroLock Knotless PEEK OPTIMA Screw-in Anchors consists of PEEK-OPTIMA (Poly-ether-ether-ketone) material. The GyroLock suture anchor comes with Dual start thread technology (fast inserting suture anchor) which allows anchor deployment in 3.5 rotation.

Vented and cannulated GyroLock anchors are designed for repair and reconstruction with soft-tissue graft. GyroLock anchor is available with multiple BioBraid suture and suture tape configurations.

GyroLock suture anchor eyelet can accommodate up to four 1.4mm tapes or four #2 sutures.

Ordering Information

Code	Product Description
6-052-01	AUXILOCK [®] GyroLock 3.5mm x 15.8mm Knotless PEEK OPTIMA Screw-In Anchor

Related Instrumentation

Code	Product Description
7-013-71	Tap for 3.5mm GyroLock Screw-In Anchor

STERILE EO