


## AUXILOCK<sup>®</sup> KNOTLESS PEEK CF PUSH-IN SUTURE ANCHOR

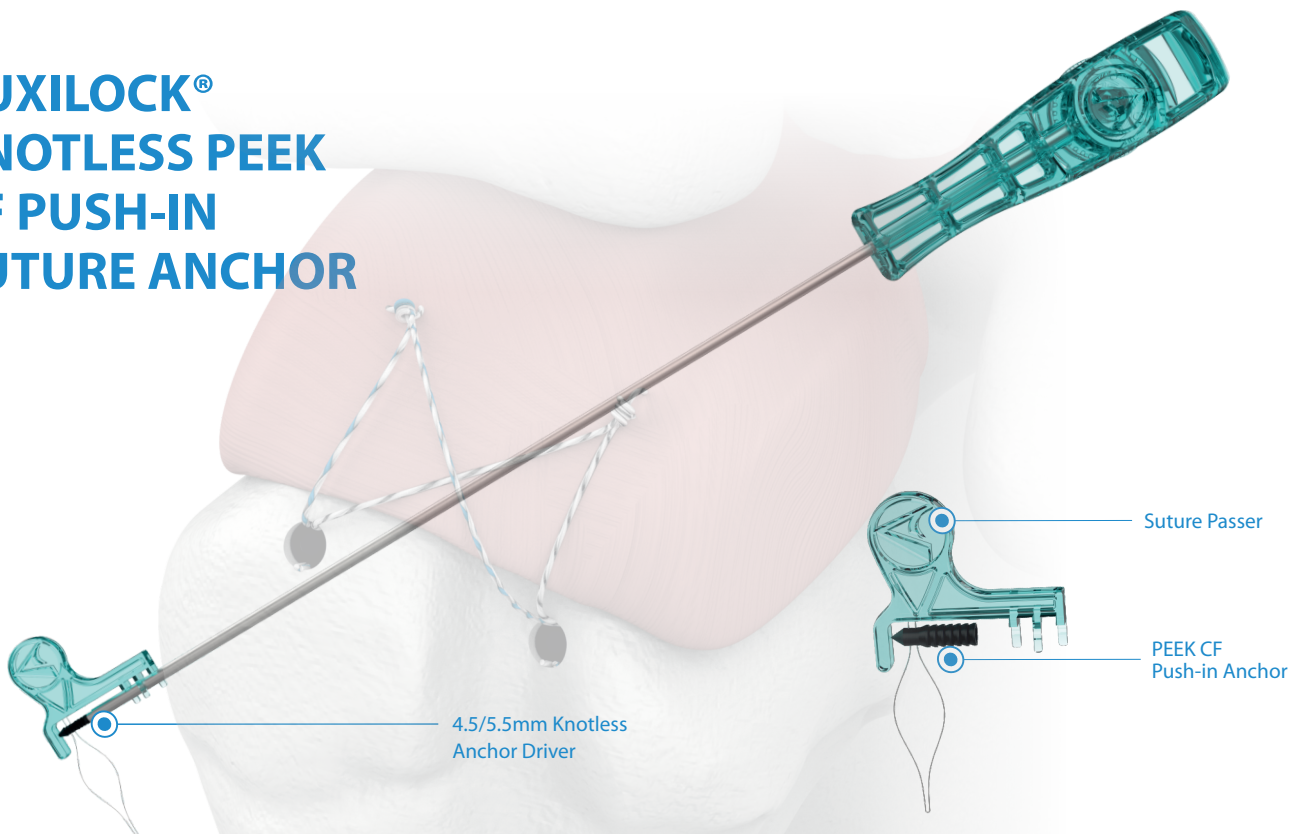

AUXILOCK<sup>®</sup> knotless push-in anchors provide a step-saving alternative to conventional knotted suture anchors. AUXILOCK<sup>®</sup> knotless push-in anchors are recommended for use in both large and small-joint soft tissue repairs. PEEK CF Push-In Suture Anchors consists of carbon fibre-reinforced PEEK OPTIMA and are available in a variety of diameters of 4.5mm and 5.5mm providing intraoperative flexibility. PEEK CF provides modulus of elasticity closely matching the cortical bone. It is non-absorbable, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation. The anchor is available with and without driver. The 4.5/5.5mm knotless anchor driver is provided separately in the instrument set. It is reusable driver that minimizes the cost of the surgery unlike the pre-loaded suture anchors.

### Ordering Information

Code	Product Description
6-003-04	AUXILOCK <sup>®</sup> 4.5mm Knotless PEEK CF Push-In Suture Anchor With Suture Passer
6-003-05	AUXILOCK <sup>®</sup> 5.5mm Knotless PEEK CF Push-In Suture Anchor With Suture Passer
6-003-14	AUXILOCK <sup>®</sup> 4.5mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer
6-003-15	AUXILOCK <sup>®</sup> 5.5mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer

### Recommended Instrumentation

Code	Product Description	Compatible Implants
7-013-23	Ø3.2mm Awl	6-003-04, 6-003-14
7-013-24	Ø4.1mm Awl	6-003-05, 6-003-15
7-013-32	4.5/5.5mm Knotless PEEK Anchor Driver	6-003-04, 6-003-05

STERILE EO